



ISQOOM®
Operating instructions





Brief overview	4
Application	6
The xCential serums	8
The xCential programmes	10
The xCential blue serums	16
The xCential blue programmes	18
Care and storage	24
Technical data and important hints	26

Thank you very much for your decision to purchase SQOOM. This high-quality concept from the medical- and dermatology expert Schick Medical offers highly effective serums and specially developed programmes for anti-ageing, whitening and cleansing, as well as applications for skin problems and back- or muscle tension.

Quite possibly the most elegant way to achieve visibly beautiful skin. The revolutionary SQOOM concept with its patented synchron technology of ultrasound and iontophoresis.

Please read these instructions carefully for optimum results and maximum treatment success.



1 | Transducer with titanium alloy – suitable for allergy sufferers

Apply the active ingredients serum generously here.

2 | Control lights

The control lights of the device show that the device is in operation (only during skin contact) by flashing regularly.

3 | On/Off POWER button

The device is switched on/off by pushing the button once.

4 | Programme selection button STEP

Any time you push the STEP button you will initiate one of the programmes.

5 | Contact handles

Always hold your SQOOM touching the contact handles on both sides.

6 | Charging station

Put the device back into the charging station after each application to keep your SQOOM ready for operation at any time.

7 | Battery control light

When the control light glows green, the device is fully charged. If it glows blue, the device is still charging.


Only use original xCential serums and original SQOOM devices for the SQOOM application. We accept no liability for results when using cosmetics or devices produced by others.



Clean your skin as usual before each application. Next, apply one of the active ingredients serums generously onto the transducer (1) or directly onto the skin area to be treated. Switch on the device (2).

Now choose the suitable application programme (3) and hold the SQOOM firmly on both sides so that you touch the contact handles (4). Then place the transducer on the skin area you want to treat and move it softly over the skin in slow circulating moves (5).

Please note: If you treat another person with your SQOOM, you should ensure that your free hand touches the other person (skin contact). This closes the circuit allowing your SQOOM to work flawlessly (blinking control light) (6).

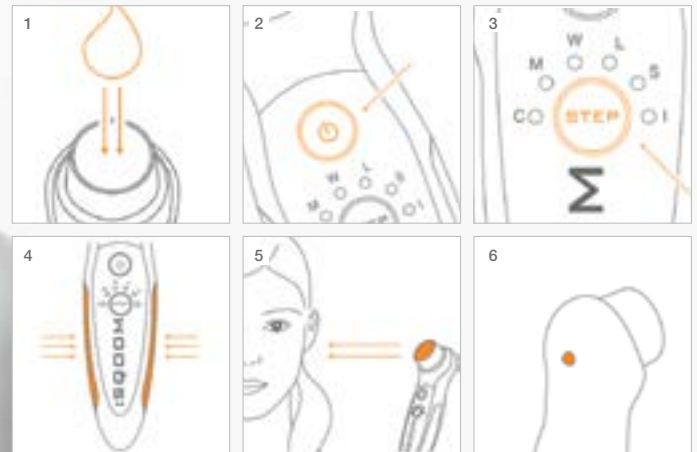

**xCential cleanser**

Interactive physical cleansing process, supplying skin with moisture and preparing it for the SQOOM treatment. Panthenol and chamomile extract can have a skin-soothing and relaxing effect.

xCential H+

The innovative H+ serum is made entirely from skin friendly natural active ingredients of natural origin for visibly beautiful and radiant skin. The main ingredient of xCential H+ is hyaluron – a renowned moisturiser. 89% of the users say “SQOOM xCential H+ gives my skin a young and beautiful appearance”*. Aloe Vera Leaf Juice gives the skin intensive moisture.

*Long term application test for 5 weeks with 75 subjects

xCential deSpot

This whitening serum contains a combination of the active ingredients sulfuraphane and genistein. This can brighten skin areas and reduce age spots in a liposomal combination.

xCential booster

24 hour concentrated moisture boost. Fresh and velvety-smooth skin with immediate effect. Wrinkle depth seems reduced after one single application and xCential booster has an 8 hour anti-wrinkle effect.

xCential acqua+

acqua+ is an enriched moisture spray, specially developed for use with all xCential serums, optimising their performance and preventing them from drying too quickly on the skin.

xCentialcleanser

- Clean your skin as usual (remove your make-up, creams and perfume).
- Apply the xCential cleanser serum generously onto the transducer or the skin area you want to treat.
- Switch on the device with the POWER button and select programme **C** with the STEP button or programme **I** (blue) for the sensitive eye area.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Wash your face with a lot of water after the cleansing treatment.
- Result: Skin cleaned down to the pores.



- We recommend carrying out the cleansing programme before each treatment.
- Apply the xCential H+ serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **W** with the STEP button or programme **I** (green) for the sensitive eye area.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- We recommend leaving the serum on your skin as a mask overnight when you treat your skin in the evening.
- If you carry out the treatment in the morning or during the day, we recommend washing off the serum after the treatment.
- Result: Visibly beautiful skin.

xCentialH+



Tip: Use programme **M** and **L** in addition to the normal H+ (**W**) application.

PROGRAMME 1 | C

Cleansing down to the pores



PROGRAMME 6 | I blue

Especially around your eyes



Treatment duration: 4 min
Application: Once a week
Use daily for quicker results.

PROGRAMME 3 | W

Reduction of wrinkles and skin firming



PROGRAMME 6 | I green

Especially around your eyes



Treatment duration: 8 - 12 min
Application: Once a week
Use daily for quicker results.

- We recommend carrying out the cleansing programme before each treatment.
- Apply the xCential deSpot serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **W** with the STEP button or programme **I** (green) for the sensitive eye area.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- We recommend leaving the serum on your skin as a mask overnight when treating your skin in the evening.
If you want to carry out the treatment in the morning or during the day, we recommend washing off the serum after the treatment.
- Result: Whitening of skin areas and age spots.



- Apply the xCential booster serum onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **S** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after the booster application and massage it into your skin with your fingers.
- Result: Gives your skin a youthful and bright look.



Tip: Use xCential acqua+ in combination before and during the application.

PROGRAMME 3 | W

Whitening application



PROGRAMME 6 | I green

Especially around your eyes



Treatment duration: 8 - 12 min
 Application: Once a week
 Use daily for quicker results.

PROGRAMME 5 | S

Concentrated moisture for 24 hours



Treatment duration: 4 min
 Application: As needed – the immediate effect for special occasions or as an ideal supplement to the xCential H+ application.



acqua+, a specially enriched moisture spray, is used as an ideal supplement to the xCential serums and xCential blue serums. It can be applied with all serums and optimises the treatment by preventing serums from drying too quickly on the skin. Additionally it supplies the skin with extra moisture with gentle active ingredients.

Spray acqua+ directly onto the skin areas you want to treat by pumping 2 – 4 times before or during the application of the xCential serums.

By using acqua+ you will use less high-value serum per treatment and the serums will last longer on your skin.





xCentialOL

Back- and muscle tension

Powerful liposomal encapsulated boswellia serrata extract is able to relax tension and relieve discomfort.

xCentialAK

Impure, greasy skin

Natural salicylic acid extracted from willow bark is able to open blocked pores. Black cohosh can have a regulating effect on the sebum production of the sebaceous glands.

xCentialRO

Reddened, sensitive skin

Magnolol and the active ingredients of caper buds can have a skin soothing and extra sensitive effect.

xCentialPS

Dry, flaky skin

Extracted from liposomal active ingredients of boswellia acid, magnolol, liposomal encapsulated liquorice active ingredients, calendula extract, dexpanthenol and pure camelia seed oil.

xCentialNA

Scars and stretch marks

Bio-active amino acids and natural onopordon acanthium extract is able to assist in the skin's own repair process.

xCentialNE

Dry, irritated skin

Magnolol, dexpanthenol, extract obtained from buds of Capparis spinosa and pure camelia seed oil can have a skin soothing and moisturising effect.

Only use original xCential serums and original SQOOM devices for the SQOOM application. We accept no liability for results when using cosmetics or devices produced by others.

xCential^{OL}

- Apply the xCential OL serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **M** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves by applying some pressure.
- Leave the remaining serum on your skin after finishing the OL application and massage it into your skin with your fingers.



xCential^{RO}

- Apply the xCential RO serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **M** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after finishing the RO application and massage it into your skin with your fingers.



PROGRAMME 2 | M

Back- and muscle tension



Treatment duration: 12 - 16 min
 Application: Daily or as needed

PROGRAMME 2 | M

Reddened, sensitive skin



Treatment duration: Face: 4 - 8 min
 Body: 12 - 16 min
 Application: 3 times a week

xCentialNA

- Apply the xCential NA serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **M** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after finishing the NA application and massage it into your skin with your fingers.



xCentialAK

- Apply the xCential AK serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **W** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after finishing the AK application and massage it into your skin with your fingers.



PROGRAMME 2 | M

Scars and stretch marks



Treatment duration: Face 4 – 8 min
 Body 12 – 16 min
 Application: Daily or as needed

PROGRAMME 3 | W

Impure, greasy skin



Treatment duration: Face 4 – 8 min
 Body 12 – 16 min
 Application: Daily or as needed

xCentialPS

- Apply the xCential PS serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **M** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after finishing the PS application and massage it into your skin with your fingers.



xCentialNE

- Apply the xCential NE serum generously onto the transducer or directly onto the skin area you want to treat.
- Switch on the device with the POWER button and select programme **M** with the STEP button.
- Hold the SQQOOM on both sides so that you touch the contact handles.
- Now put the transducer onto the skin area you want to treat and move it on your skin with slow, soft circulating moves.
- Leave the remaining serum on your skin after finishing the NE application and massage it into your skin with your fingers.



PROGRAMME 2 | M

Dry, flaky skin



Treatment duration: Face 4 – 8 min
Body 8 – 16 min
Application: Daily or as needed

PROGRAMME 2 | M

Dry, irritated skin



Treatment duration: Face 4 – 8 min
Body 12 – 16 min
Application: Daily or as needed

Please read the following tips to help you get maximum enjoyment from your device.

Only use original xCential serums and original SQOOM devices for the SQOOM application. We accept no liability for results when using cosmetics or devices produced by others.



Extended international warranty

Please register your device for the FREE extended 2-year complete warranty and become SQOOM VIP. Extending the legally prescribed warranty and liability, you can use the extended worldwide complete SQOOM warranty free of charge and receive many other benefits. www.sqoom.com/en/meta/vip-registration/

Customer service number: +49 (0) 741 174 124 0
Monday to Friday from 8am - 6pm

- Store in a dry place and do not immerse in water.
- Clean the transducer with a moist cloth.
- Put the device in the charging station after switching it off and make sure that the blue LED light glows.
- The blue LED light turns to green when the charging is finished. If there are problems with charging, the LED will blink blue.
- If the battery of the device is discharged completely, please charge it for at least three hours before using the device again.
- When completely charged, the device has approximately four hours of usage time. One full charge can provide over two weeks of use at 12 minutes per day.
- Press the STEP button for about three seconds to display the remaining battery charge. If the LEDs of the 6 programmes glow, the device will be charged completely. If only two LEDs glow, there is approximately 30 minutes of charge left. When the signal tone sounds, the battery must be charged.
- If you do not use the device within two minutes after switching it on, it will switch off automatically.
- The programmed application time for each step is four minutes. When the four minutes have expired, a signal tone sounds. After this, the device switches to stand-by mode. If it is not switched on again within two minutes, it will switch off automatically.
- The rated input voltage is 100 V to 240 V AC. So the device can also be taken to your journeys and used worldwide.
- Remove the power plug of the charging station as soon as your device is fully charged.
- Keep all electric contacts clean, dust-free and dry.

Device Designation

Device SQOOM m2

Charging Station

Mains Voltage AC 100 - 240V, 50/60 Hz

Output/ Adapter DC 4.2 V, 1A

Power 7 W

Weight 133 g

Hand Device

Lithium/ ion battery

Output/ Adapter DC 3,7V / 2200mA

Weight 194 g (gross)

Ultrasound

Frequency 1 MHz

Intensity levels from 0.1 W/ cm² to 0.5 W/ cm²

Ion

Frequency 1,150Hz ~ 1,785Hz



In the interest of recycling raw materials, the device should not be disposed of with the normal household waste at the end of its operating life. The disposal can be executed via Schick Medical GmbH or at applicable collecting points which you can find out from your local council or municipal administration.

Only use original xCential serums and original SQOOM devices for the SQOOM application. We accept no liability for results when using cosmetics or devices produced by others.

Please do not treat one skin area longer than 20 minutes per day.

Do not treat an area with SQOOM which is affected by skin cancer or surgery sites.

Do not use SQOOM if you have a cardiac pacemaker.

In cases of hyperthyroidism, do not treat the throat or chest area with SQOOM.

Do not treat the area around joint implants with SQOOM. This might result in increasing pain.

Do not treat the abdomen or chest area with SQOOM during pregnancy

Avoid the contact of serum and transducer with your eyeballs, eye lids and mucous membranes.

Keep the device and serum/cream out of the reach of children.

Do not store the serum/cream at more than 25°C.

The serums/creams are exclusively made for external skin application.

If you treat another person with your SQOOM, you should ensure that your free hand touches the other person (skin contact). This will close the circuit and ensure the proper operation of your SQOOM.

© All rights are reserved to Schick Medical GmbH. Any copies or editing beyond the copyright framework are not permitted without prior written admission.

Schick Medical GmbH,
D-78628 Rottweil



A registered trademark of Schick Medical GmbH

Auf dem Wall 29

D-78628 Rottweil

Tel.: +49(0)741/174 124 0

Fax: +49(0)741/174 124 19

info@SQOOM.com

www.SQOOM.com